

staan®  
— MADE FOR THE WORLD —



ISO 13485 : 2016



Inspire  
*options unlimited.....*

C-Arm Compatible  
Multipurpose Electrical Table



Multipurpose electro mechanical operating table is ergonomically designed to bring out the versatility in the Operation theatre as per the surgical requirement. It is mobile, electrical operated surgical table for Ortho, Neuro, Spine, Cardio thoracic, Gastro Enterology, Bariatric, PCNL, Urology, Gynec, Lapro, Ophthalmic, ENT and General Surgeries



Electrical



Function



Zero function



Electrical Brakes



Override Panel



Battery Backup



Incorporated kidney elevator



C-Arm compatible top leg traction attachment



C-Arm compatible URS, TURP and Gynecological surgical positions



Cardiothoracic surgical position

- Modular and stylish design
- Smooth and jerk free movements
- Versatile Usage Providing Maximum Comfort To The Patient And Surgeon
- Blunted corners and smooth edges assures patient and user safety.
- Side rails on both sides
- 250Kgs Dynamic weight bearing capacity
- All the material are made of SS 304 and base is covered by molded ABS cover for aesthetic look
- Scratch free Imported poly urethane twin disc castors
- Self lubricating jerkless, weldless jointless steel column
- Sealed bottom prevents water and dust from entering inside
- The table operates on 220V/5 Amps domestic power supply
- Controllers, handsets, cables, actuators, override panels and batteries are all world class



Shoulder arthroscopy surgical position



Pateint in sitting, neuro surgical position



C-Arm compatible knee-chest spine position



Lapro and Gynae Surgeries




Technical Data

- Option 1
- Option 2
- Option 3

	OPTION 1	OPTION 2	OPTION 3
Mode	Normal height Electrically operated	Low height Electrically operated	Ultra low height Electrically operated
Specialization	<b>MULTIPURPOSE</b> Compatible for all surgeries including orthopaedic attachment	<b>MULTIPURPOSE</b> Compatible for all surgeries excluding orthopaedic traction attachment (due to height constraint)	<b>LAPRO-GYNEC</b> compatible for all Lapro Gynec surgeries

Electrical Movements

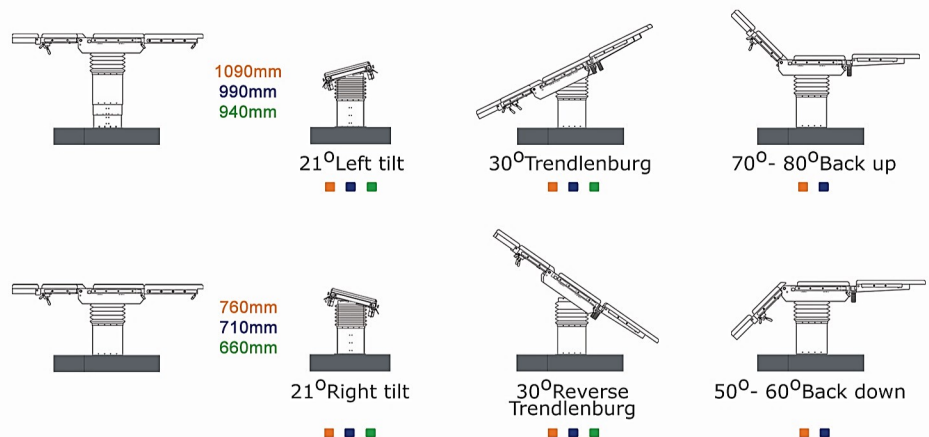

Manual Movements


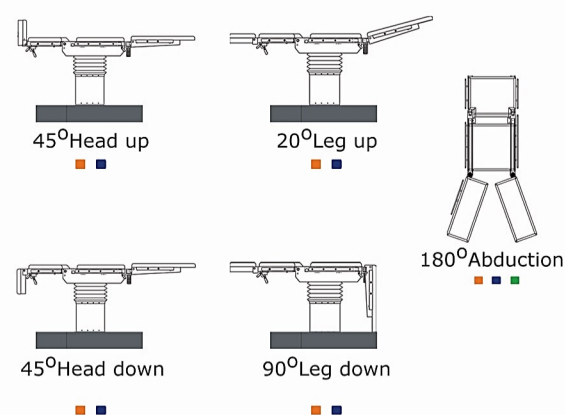


Table Top Dimensions

Length  
2000mm

Width  
550mm without side rails  
600mm with side rails

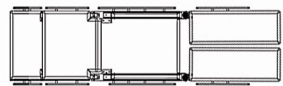

Watch demo video @  
"Staan" channel

